



Endo Agar

M029

Intended Use

Endo Agar is a selective medium recommended for confirmation of the presumptive test for members of the coliform group from clinical and non-clinical samples.

Composition**

Ingredients	Gms / Litre
Peptone	10.000
Lactose	10.000
Dipotassium phosphate	3.500
Sodium sulphite	2.500
Basic fuchsin	0.500
Agar	15.000
Final pH (at 25°C)	7.5±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 41.5 grams in 1000 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Mix well before pouring into sterile Petri plates. If the solidified culture medium is somewhat too red, then to remove the colour add a few drops (max. 1 ml/litre) of a freshly prepared 10% Sodium sulphite solution and boil.

Principle And Interpretation

Endo Agar was developed by Endo to differentiate gram-negative bacteria on the basis of lactose fermentation, while inhibiting gram-positive bacteria (2). Inhibition of the later was achieved without the use of bile salts as was traditionally used. Endo was successful in inhibiting gram-positive bacteria on his medium by the incorporation of sodium sulphite and basic fuchsin. The resulting Endo Agar, also known as Fuchsin Sulphite and Infusion Agar, was used to isolate the typhoid bacilli. Many modifications of this media have been done over the years. Endo Agar is recommended by APHA as an important medium in the microbiological examination of water and wastewater, dairy products and foods (1,5,6). Endo Agar is used to confirm the detection and enumeration of coliform bacteria following presumptive test of drinking water. It is also used for the detection and isolation of coliforms and faecal coliforms from milk, dairy products and food.

The medium contains peptone which provide nitrogen, carbon, vitamins and minerals required for bacterial growth. Sodium sulphite and basic fuchsin make this medium selective by suppressing gram-positive organisms. Coliforms produce pink colonies on fermentation of lactose while lactose non-fermenters produce colourless colonies on the medium.

With *Escherichia coli*, this reaction is very pronounced as the fuchsin crystallizes, exhibiting a permanent greenish metallic luster (fuchsin luster) to the colonies. Medium should be stored away from light to avoid photo-oxidation.

Type of specimen

Clinical samples - urine,; Food and dairy samples; Water samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (3,4).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (5,6).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards.(1)

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

Limitations :

1. Besides *Enterobacteriaceae*, other gram negative bacteria and yeasts may also grow.
2. Avoid exposure of the medium to light, as it may lead to photo oxidation and decrease productivity of the medium.
3. Overheating of the medium must be avoided, as it may destroy the productivity of the medium.
4. Further biochemical tests must be carried out for further confirmation.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light pink to purple homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Orangish pink coloured, clear to slightly opalescent gel with fine precipitate forms in Petri plates.

Reaction

Reaction of 4.15% w/v aqueous solution at 25°C. pH : 7.5±0.2

pH

7.30-7.70

Cultural Response

M029: Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of Colony
Cultural Response				
<i>Bacillus subtilis subsp. spizizenii</i> ATCC 6633 (00003*)	≥10 ⁴	inhibited	0%	
* <i>Klebsiella aerogenes</i> ATCC 29212 (00087*)	50-100	good-luxuriant	≥50%	pink
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50-100	none-poor	≤10%	pink, small
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant	≥50%	pink to rose red with metallic sheen
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	50-100	good-luxuriant	≥50%	pink, mucoid
<i>Proteus vulgaris</i> ATCC 13315	50-100	good-luxuriant	≥50%	colourless to pale pink
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50-100	good-luxuriant	≥50%	colourless, irregular
<i>Salmonella Typhi</i> ATCC 6539	50-100	good-luxuriant	≥50%	colourless to pale pink
<i>Staphylococcus aureus subsp. aureus</i> ATCC 25923 (00034*)	≥10 ⁴	inhibited	0%	
<i>Enterobacter cloacae</i> ATCC 13047 (00083*)	50-100	good	40-50%	pink
<i>Salmonella Typhimurium</i> ATCC 14028 (00031*)	50-100	good-luxuriant	≥50%	colourless
<i>Salmonella Enteritidis</i> ATCC 13076 (00030*)	50-100	good-luxuriant	≥50%	colourless
<i>Shigella flexneri</i> ATCC 12022 (00126*)	50-100	good-luxuriant	≥50%	colourless

#- Formerly known as *Enterobacter aerogenes*

Key : *Corresponding WDCM numbers.

Please refer disclaimer Overleaf.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Use before expiry date on the label.

Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (3,4).

Reference

1. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
2. Endo S., 1904, Zentralbl. Bakteriol., Abt. 1, Orig.35:109-11
3. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
4. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
5. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.


Revision : 04/ 2019

 IVD

In vitro diagnostic medical device



CE Marking

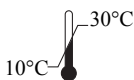

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U, Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.