



Millon's Reagent

R097

Intended Use:

Millon's reagent is an analytical reagent used to detect the presence of tyrosine, the only amino acid containing a phenol group.

Composition**

Ingredients

Mercuric Nitrate	160.0g
Mercurous Nitrate	160.0g
Conc. Nitric acid	400.0ml
Distilled water	600.0ml

**Formula adjusted, standardized to suit performance parameters

Directions

- 1) Add 1ml of Millon's reagent to 1-2ml of test solution.
- 2) If desired heat the solution till it boil.
- 3) Observe the pink to brick red precipitate or solution which indicate positive test.

Principle And Interpretation

In Millon's test, the phenol group of tyrosine is first nitrated by nitric acid in the test solution. Then the nitrated tyrosine complexes mercury (I) and mercury (II) ions in the solution to form a pink to brick red precipitate or a solution. Some proteins containing tyrosine initially forms a white precipitate that turns red when heated, while others form a red solution immediately. Both results are considered positive.

Warning and Precautions

In Vitro diagnostic use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations

1) Millon's test is specific for phenolic compounds, so that other proteins must be confirmed by other tests for proteins such as the biuret test and the ninhydrin reaction.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Colourless to light greenish yellow solution

Clarity

Clear solution with no insoluble particles

Test

Add 1ml of Millon's reagent to 1-2ml of test solution. If desired heat the solution till it boil.

Results

Pink to brick red precipitate or solution	presence of proteins (tyrosine, the only amino acid containing a phenol group)
No colour change	absence of protein

Storage and Shelf Life

Store between 10-30°C in tightly closed container and away from bright light. Use before expiry date on label. On opening, product should be properly stored in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (2, 3).

Reference

1)Walsh, Edward O'Farrell (1961). An Introduction to Biochemistry. London: The English Universities Press Ltd. pp. 406–407. OCLC 421450365.

3) Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.

4) Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock, D.W. (2015) Manual


Revision : 00/ 2017



In vitro diagnostic medical device



CE Marking

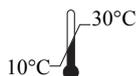

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
B /4-6 , MIDC, Dindori, Nashik MH

www.himedialabs.com



CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,

www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.