

Complete Supplement Mixture w/o LYS

G108

Complete Supplement Mixture w/o LYS is used as a dropout supplement for all strains of *Saccharomyces cerevisiae*.

Composition** :

Ingredients	Milligrams/Litre
Adenine	10.00
L-Arginine HCl	50.00
L-Aspartic acid	80.00
L-Histidine HCl	20.00
L-Isoleucine	50.00
L-Leucine	100.00
L-Methionine	20.00
L-Phenylalanine	50.00
L-Threonine	100.00
L-Tryptophan	50.00
L-Tyrosine	50.00
Uracil	20.00
L-Valine	140.00

** Formula adjusted, standardized to suit performance parameters

Directions :

Suspend 740 mg in 1000 ml distilled water. Sterilize by autoclaving at 10 lbs pressure (115°C) for 20 minutes. Mix well and dispense as desired.

Principle and Interpretation :

Complete Supplement Mixture w/o LYS is used as a dropout supplement for all strains of *Saccharomyces cerevisiae*. Complete Supplement Mixture (CSM) is composed of all the amino acids required for the vigorous growth of *Saccharomyces cerevisiae*. It also contains uracil. CSM w/o LYS does not contain lysine. This makes it a dropout growth medium for yeast cells. It can be mixed with yeast nitrogen base (YNB), ammonium sulphate and an appropriate carbon source (glucose/galactose/raffinose) to produce a media suitable for the growth of lysine prototrophs and wild type strains of *S. cerevisiae* but lysine auxotroph cannot grow in this media. Therefore, CSM w/o LYS can be used as drop-out formulation for the selection of auxotrophic requirements and transformants.

Please refer disclaimer Overleaf

1



Registered Office :

23, Vadhani Industrial Estate, LBS Marg,
Mumbai - 400 086, India.
Tel. : (022) 4017 9797 / 2500 1607
Fax : (022) 2500 2286

Commercial Office

A-516, Swastik Disha Business Park,
Via Vadhani Indl. Est., LBS Marg,
Mumbai - 400 086, India

Tel: 00-91-22-6147 1919
Fax: 6147 1920, 2500 5764
Email : info@himedialabs.com
Web : www.himedialabs.com

Quality Control :

Appearance of Powder :

White to off-white, homogeneous, free flowing powder.

Colour and Clarity :

Colourless, clear solution without any precipitate.

Cultural Response :

Cultural characteristics observed after an incubation at 25-30°C for 18 - 48 hours.

Organisms (ATCC)


Saccharomyces cerevisiae ATCC 9763

Growth

good-luxuriant

Storage and Shelf-life :

Store below 30°C and the prepared medium at 2 - 8°C. Use before expiry date on the label.

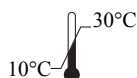

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India

PIG108_1/0619

G108-01

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal diagnostic or therapeutic use but for laboratory, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.

HiMedia Laboratories Pvt. Ltd. A-516, Swastik Disha Business Park, Via Vadhani Ind. Est., LBS Marg, Mumbai-400086, India.
Customer care No.: 022-6147 1919 Email: techhelp@himedialabs.com Website: www.himedialabs.com