



M-CP Agar Base

M1354

Intended Use:

Recommended by the directive council of the European Union 98/83/EC for isolation and enumeration of *Clostridium perfringens* from water samples using membrane filtration technique.

Composition**

Ingredients	Gms / Litre
Tryptose	30.000
Yeast extract	20.000
Sucrose	5.000
L-Cysteine hydrochloride	1.000
Magnesium sulphate, 7H ₂ O	0.100
Bromocresol purple	0.040
Ferric chloride, 6H ₂ O	0.090
Indoxyl β-D-glucoside	0.060
Agar	15.000
Final pH (at 25°C)	7.6±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 35.60 grams (the equivalent weight of dehydrated medium in 485 ml) in 485 ml purified/distilled water. Heat to boiling to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50° C. Aseptically add the rehydrated contents of 1 vial of M-CP Selective Supplement I (FD153) and 1 vial of M-CP Selective Supplement II (FD154) or rehydrated contents of 1 vial of M-CP Selective Supplement II, Modified (FD154A). Mix well and pour into sterile Petri plates.

Principle And Interpretation

Clostridial species are one of the major causes of food poisoning/ gastro-intestinal illnesses. They are gram-positive, spore-forming rods that occur naturally in the soil (3). Among the family are: *Clostridium botulinum* which produces one of the most potent toxins in existence; *Clostridium tetani*, causative agent of tetanus; and *Clostridium perfringens* commonly found in wound infections and diarrhoea cases. The use of toxins to damage the host is a method deployed by many bacterial pathogens. The major virulence factor of *C. perfringens* is the CPE enterotoxin, which is secreted upon invasion of the host gut, and contributes to food poisoning and other gastrointestinal illnesses (3). Several solid media have been devised for quantitation of *C. perfringens*. The selectivity of the media is achieved by incorporation of one or more antibiotics that inhibit certain anaerobes or facultative anaerobes.

M-CP Agar Base is prepared as per the formula of Armon and Payment (1). It is also recommended by the Directive of the Council of the European Union 98/83/EC (4) for isolation and enumeration of *Clostridium perfringens* from water sample using membrane filtration technique.

Tryptose, yeast extract provide nitrogenous, carbonaceous compounds, long chain amino acids, vitamins, minerals and other essential growth factors while sucrose is the fermentable carbohydrate. Bromocresol purple serves as a pH indicator. Indoxyl-β-D-glucoside is a chromogenic substrate for β-D-glucosidase or cellobiase and phenolphthalein diphosphate for the detection of acid phosphatase. The addition of D-cycloserine and polymyxin B (FD153) makes the medium inhibitory to accompanying non-clostridial microflora and thus allows analysis of both clostridial vegetative cells and spores. Further selectivity is provided by incubation under anaerobic conditions. Yellow (cellobiase-negative) colonies becoming old rose to pink-red upon exposure to ammonia fumes for 30 seconds are considered to be presumptive *C. perfringens*. Colour differentiation on M-CP Agar Base is sometimes difficult, so typical colonies (yellow turning into pink) as well as atypical

colonies (green or those that remain yellow upon exposure to ammonia fumes) are picked for confirmation. Presumptive *C. perfringens* can be confirmed by sulphite reduction, gram-positive, sporulating rods, non-motile, reduction of nitrate, gelatine liquefaction, lactose fermentation and other biochemical tests (8).

Type of specimen

Clinical samples; Water samples; Food samples

Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (2).

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (7).

After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions

In Vitro diagnostic Use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets

Limitations

1. Colour differentiation on M-CP Agar Base is sometimes difficult, so typical colonies (yellow turning into pink) as well as a typical colonies (green or those that remain upon exposure to ammonia fumes) are picked for confirmation.
2. Due to variable nutritional requirements, some strains may show poor growth on this medium.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Light yellow to light green homogeneous free flowing powder

Gelling

Firm, comparable with 1.5% Agar gel

Colour and Clarity of prepared medium

Purple coloured clear to slightly opalescent gel forms in Petri plates

Reaction

Reaction of 7.12% w/v aqueous solution at 25°C. pH : 7.6±0.2

pH

7.40-7.80

Cultural Response

Cultural characteristics observed after an incubation at 44°C for 24-48 hours with added contents of 1 vial of M-CP Selective Supplement I (FD153) and 1 vial of M-CP Selective Supplement II (FD154) or rehydrated contents of 1 vial of M-CP Selective Supplement II, Modified (FD154A) under anaerobic conditions.

Organism	Growth	Recovery	Inoculum (CFU)	Phosphatase test(on exposure to ammonia
<i>Clostridium perfringens</i> ATCC 12924	good-luxuriant	>=50%	50-100	positive, yellow coloured colonies(colonies become old rose to light pink-red on exposure to ammonia fumes for 30 seconds)

<i>Clostridium perfringens</i> ATCC 13124 (00007*)	good-luxuriant	$\geq 50\%$	50-100	positive, yellow coloured colonies(colonies become old rose to light pink-red on exposure to ammonia fumes for 30 seconds)
<i>Clostridium perfringens</i> ATCC 10543 (00174*)	good-luxuriant	$\geq 50\%$	50-100	positive, yellow coloured colonies(colonies become old rose to light pink-red on exposure to ammonia fumes for 30 seconds)
<i>Clostridium perfringens</i> ATCC 12916 (00080*)	good-luxuriant	$\geq 50\%$	50-100	positive, yellow coloured colonies(colonies become old rose to light pink-red on exposure to ammonia fumes for 30 seconds)
<i>Clostridium bifermentans</i> NCTC 506 (00079*)	good-luxuriant	$\geq 50\%$	50-100	negative, blue coloured colonies
<i>Escherichia coli</i> ATCC 8739 (00012*)	inhibited	0%	$\geq 10^3$	
<i>Escherichia coli</i> ATCC 25922 (00013*)	inhibited	0%	$\geq 10^3$	
<i>Salmonella</i> Typhi ATCC 6539	inhibited	0%	$\geq 10^3$	
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	inhibited	0%	$\geq 10^3$	

Key ; (*) Corresponding WDCM numbers

Storage and Shelf Life

Store between 15-25°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

Reference

1. Armon R. and Payment P., 1988, Can. J. Microbiol., 34:78-79.
2. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
3. Czeczulin J. R., Hanna P. C., McClane B. A., 1993, Infect. Immun., 61: 3429-3439.
4. Directive of the Council of the European Union 98/83/EC
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.

7. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2001, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
8. Sartory D. P., Field M., Curbishley S. M., Pritchard A. M., 1998, Lett. Appl. Microbiol., 27:323-327.


Revision : 04 / 2018



In vitro diagnostic medical device



CE Marking

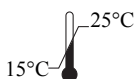

Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,
23 Vadhani Industrial Estate,
LBS Marg, Mumbai-86, MS, India



CE Partner 4U ,Esdoornlaan 13, 3951
DB Maarn The Netherlands,
www.cepartner4u.eu

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.