



## HiColiform™ Broth, Modified

M1850

### Intended use

Recommended for the detection and confirmation of *Escherichia coli* and total coliforms from water samples, using a combination of chromogenic and fluorogenic substrates. It can also be used for clinical samples.

### Composition\*\*

Ingredients	Gms / Litre
Peptone	5.000
Sodium chloride	5.000
Potassium sulfate	1.000
Dipotassium hydrogen phosphate	4.000
Potassium dihydrogen phosphate	1.000
Sodium lauryl sulphate	0.100
Sodium puruvate	1.000
Chromogenic substrate	0.100
Fluorogenic substrate	0.100
IPTG	0.100
Final pH ( at 25°C)	6.8±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 17.4 grams in 1000 ml purified/distilled water. Heat if necessary to dissolve the medium completely. Mix well and dispense in tubes or flasks as desired. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to 45-50°C.

### Principle And Interpretation

HiColiform™ Broth, Modified was designed for detection and confirmation of *Escherichia coli* and total coliforms from water samples using a combination of chromogenic and fluorogenic substrates. *Escherichia coli* can be distinguished from other coliforms by its unique ability to fluoresce in the presence of fluorogenic substrate (2, 3). The fluorogenic substrate is split by enzyme  $\beta$ -glucuronidase especially present in *Escherichia coli*. The reaction is indicated by the development of a blue fluorescence under UV light. The presence of total coliforms is indicated by blue-green colourations due to the cleavage of the chromogenic substrate. IPTG amplifies enzyme synthesis and increases the activity of  $\beta$ -galactosidase. Peptone provides carbon, nitrogen substances, long chain amino acids, vitamins and other essential growth nutrients and is useful for the simultaneous detection of indole production. The phosphate salts provide buffering action for rapid growth of coliforms. Sodium chloride helps to maintain the osmotic balance. Sodium lauryl sulphate makes the medium selective by inhibiting accompanying microflora, especially the gram-positive organisms.

### Type of specimen

Clinical samples - urine; Water samples.

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (1).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions :

In Vitro diagnostic Use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

## Limitations

1.  $\beta$ -glucuronidase is present in 97% of *E.coli* strains, however few *E.coli* may be negative.
2. Some species may show poor growth due to nutritional variations.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Cream to yellow homogeneous free flowing powder

### Colour and Clarity of prepared medium

Light yellow coloured, clear to slightly opalescent solution in tubes

### Reaction

Reaction of 1.74% w/v aqueous solution at 25°C. pH : 6.8±0.2

### pH

6.60-7.00

### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Colour of medium	Fluorescence (under uv)
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50-100	luxuriant	blue-green	negative reaction
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	luxuriant	blue-green	positive reaction

Key : (#) Formerly known as *Enterobacter aerogenes*. (\*) Corresponding WDCM numbers

## Storage and Shelf Life

Store between 15-25°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

## Reference

1. Baird R.B., Eaton A.D., and Rice E.W., (Eds.), 2015, Standard Methods for the Examination of Water and Wastewater, 23rd ed., APHA, Washington, D.C.
2. Feng P.C.S. and Hartman P.A., 1982, J. Appl. Environmental Microbiol. 43. 1320-1323.
3. Harsen W., and Yourassowsky, 1984, J. Clin. Microbiol. 20. 1177-1179.
4. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision :03 / 2020



In vitro diagnostic medical device



CE Marking

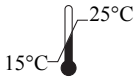


Storage temperature



Do not use if package is damaged



HiMedia Laboratories Pvt. Limited,  
23 Vadhani Industrial Estate,  
LBS Marg, Mumbai-86, MS, India

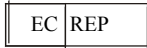

CE Partner 4U ,Esdoornlaan 13, 3951  
DB Maarn The Netherlands,  
[www.cepartner4u.eu](http://www.cepartner4u.eu)

**Disclaimer :**

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.